

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 1 OF 43

STRmix™ v 2.7 Probabilistic Genotyping Software Operating Instructions

1 Guiding Principles and Scope

- 1.1 This procedure describes the use of STRmix™ v2.7 for the interpretation of PowerPlex® Fusion DNA profiles run on 3500xL or 3130xL Genetic Analyzers within the NYC OCME Department of Forensic Biology. Readers are also referred to the STRmix™ v2.7 Users and Operation manuals for additional information.
- 1.2 For STRmix™ set-up instructions please refer to [QC702a – STRmix™ v2.7 Set-Up Instructions](#).
- 1.3 For STRmix™ set-up instructions for the QC Monitoring Program refer to [QC702b](#).

2 Preparing Data for a STRmix™ Analysis

- 2.1 Before performing your STRmix™ analysis, the following actions must be taken:
 - 2.1.1 Verify that the sample is suitable for STRmix™ analysis (refer to [Interpretation of PowerPlex® Fusion data run on 3500xL](#)).
 - 2.1.2 Evaluate your replicates, if applicable. If there are drastic inconsistencies with the alleles present between replicates and/or one has little information, only the amplification with the most information should be used, or a third amplification may be warranted.
 - 2.1.3 Determine the best described Number of Contributors to the sample (NOC). Refer to [Interpretation of PowerPlex® Fusion data run on 3500xL](#) regarding the procedure for determining the number of contributors.
 - 2.1.4 Create folder(s) with the FB, FBS, or F case number for the STRmix™ runs within the STRmix Data folder: M:\STR_Data\STRmix Data.
 - e.g. FB16-01234, FBS16-05678 or F25-01234
 - 2.1.5 Confirm that the STR data is prepared correctly for STRmix™ analysis:
 - 2.1.5.1 Evidence samples must only be amplified in PowerPlex® Fusion in order to undergo STRmix™ analysis. This procedure is specifically for evidence amplified in PowerPlex® Fusion data run on the 3500xL or 3130xL Genetic Analyzers at the NYC OCME Department of Forensic Biology.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 2 OF 43

- 2.1.5.2 Sample data must be assembled into the appropriate format for STRmix™ input. The standard input for STRmix™ is a .txt file. See the [GeneMarker v3.0 Operation Manual](#) for instructions on exporting data for STRmix™ input.
- 2.1.5.3 Evidence samples must be edited to remove all artifacts, including pull-ups, spikes, dye blobs, etc. before being input into STRmix™. **Back and forward stutter in 3130xL data** should not be removed, and **back, forward, half back, and double back stutter in 3500xL data** should not be removed before importing into STRmix™. Refer to [STR Results Interpretation - PowerPlex® Fusion & STRmix™](#), [Interpretation of PowerPlex® Fusion data run on 3500xL](#), and the [Appendix for PowerPlex® Fusion Stutter](#).
- 2.1.5.4 Reference samples must be edited to remove all artifacts including stutter. Incomplete or tri-allelic loci must not be imported into STRmix™ for a reference sample – remove all allele(s) for that locus within the text file. If a possible drop-in peak is present in a reference sample, remove this peak from the text file before STRmix™ import.
- 2.1.5.5 Non-numeric values such as OL or OB, < or > are not permitted within the STRmix™ input files. Unambiguous alleles including those that are rare should appear in the corresponding input file as their actual allelic size designation, for example D21: 30.1. If an actual allelic size designation cannot be determined, the data for this locus should be removed completely from the text file and the locus should be ignored.
- 2.1.5.6 To modify a STRmix™ input text file: open the STRmix™ .txt file associated with the appropriate STR project (e.g. in Notepad or Microsoft Excel®). Locate the sample and locus containing the non-numeric value within the .txt file and manually replace the value with the appropriate actual allelic size designation. Save the .txt file **replacing the original file**.
- 2.1.5.7 An attempt should be made to amplify reference samples in PowerPlex® Fusion (see [Case Management](#)). If unavailable, STRmix™ allows the user to calculate a likelihood ratio when the evidence and reference samples are analyzed in different autosomal typing kits. LR's will only be calculated for those loci in common between the two kits.
- 2.1.5.7.1 If a reference sample was not amplified in Fusion, the data must be converted to a .txt file for import into STRmix™ with the locus order matching that of the evidence (PowerPlex® Fusion order). The reference sample data can be converted to the proper PowerPlex® Fusion order and appropriate .txt file format using the following macro: [Identifiler to Fusion Exemplar STRmix Input Creation](#).
- 2.1.5.8 If a DNA donor is being used from one sample to condition or compare to another, use the following macro: [Reference profile for STRmix Input Creation](#).

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 3 OF 43


2.1.6 Evaluate evidence samples to determine if a locus needs to be ignored before STRmix™ deconvolution is performed. A comment should be added to the **Case Notes** field within the STRmix™ analysis to indicate why the locus was ignored for that run.

2.1.6.1 The following is a list of reasons data may need to be ignored at a particular locus. For any situation not covered in the list below, the technical leader should be consulted.

- Tri-allelic pattern
- Unresolved allelic or stutter peak that is visible above the AT
- OB/OL allele or stutter peak that cannot be assigned a correct allelic designation
- Stutter or allelic peak for an allele belonging to locus A is being called in a neighboring locus B (ignore both loci)
- Where a conditioning sample does not have data at a locus that is present in the evidence sample

2.1.7 Loci should not be ignored for likelihood ratio calculations due to a partial comparison sample. Ensure that your reference sample text file has been updated appropriately as described above in 2.1.5

2.2 Launch the STRmix™ application and prepare the scenarios to be run in STRmix™. Open the STRmix™ software by locating STRmix™ in the task bar or by double clicking on the STRmix™ icon on the desktop. The main menu is shown below:

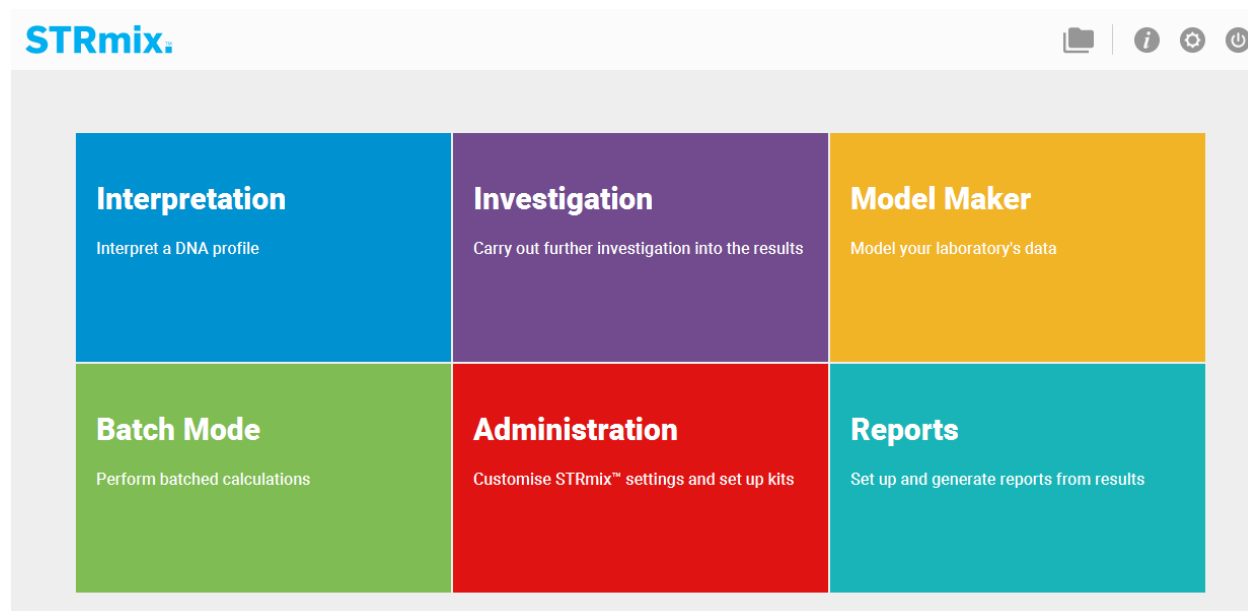

2.2.1 For deconvolutions of evidence profiles (Interpretation), go to Section 3.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS


STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 4 OF 43

2.2.2 For comparisons to previously deconvoluted profiles (Investigation: LR From Previous), go to Section 4.

2.2.3 To set up multiple STRmix™ runs to run consecutively in Batch Mode, go to Section 6.

3 Deconvolutions with STRmix™ (Interpretation)

3.1 Select **Interpretation** from the **Main Menu**. This will open the Interpretation Setup screen:

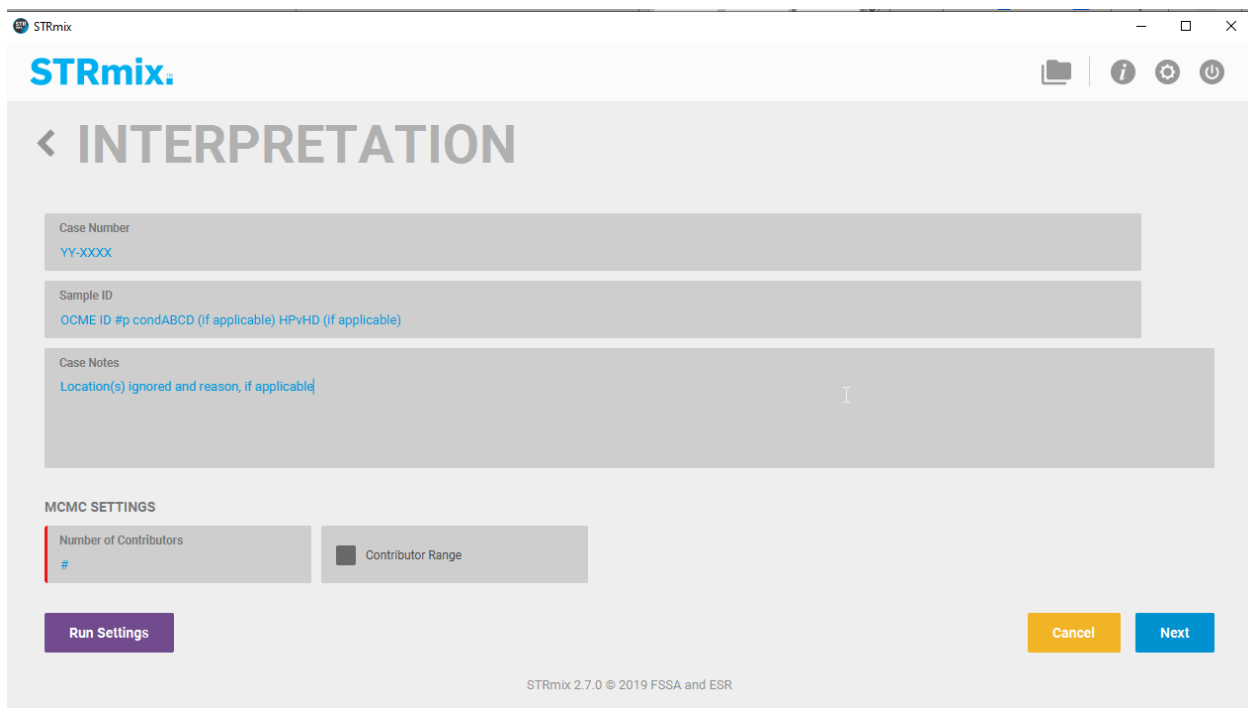

3.2 The STRmix™ output folder and file names are created by stringing together the values in the **Case Number** and **Sample ID** fields in the software followed by the date and time of the analysis run. The information in the **file name** is separated by dashes. Therefore, if other characters are entered, such as a comma, underscore, period, etc., the software will convert them into dashes.

3.3 Refer to [Interpretation of PowerPlex® Fusion data run on 3500xL](#) for guidance on when a **conditioned contributor** may be applied. A deconvolution of the evidence sample without conditioning and an LR against a potential conditioned contributor may need to be run first.

3.4 An LR may be run in conjunction with the deconvolution in the following scenarios:

- Single source evidence sample that did not require a STRmix™ deconvolution for determination of a profile, needing an LR to a matching comparison sample.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 5 OF 43

- Evidence mixture sample undergoing deconvolution and an LR check to determine if a reasonably expected reference sample can be used for further conditioning (ex. car owner on a swab from the steering wheel).
- To assess a probative comparison sample's presence within a mixture where that reference sample data is contained within the same evidence file (ex. victim's blood on suspect's clothing.)

3.5 The following naming conventions should be used for deconvolutions:

3.5.1 Evidence Files

- **Case Number** = YY-XXXXX (do not include "FB" or "F")
- **Sample ID** = item and sample identifiers consistent with evidence electropherogram. Suffixes such as '_mcon' or '_reamp' should not be included. Exclude additional identifiers that are specific to sample reps or recuts.
- **NOC**: #p
- **Conditioned**: condInitials (if applicable)
- **LR**: scenario for LR (if applicable, refer to 3.5.3)
- **Case notes** = a comment should be added here if a locus is ignored, indicating the reason: e.g. "D2S441 was ignored due to an unresolved allelic peak"

3.5.2 Suspect Files

- **Case Number** = SYY-XXXXX (do not include "FB" or "F")
- **Sample ID** = item and sample identifiers consistent with evidence electropherogram. (include evidence file FB# without the "FB") Suffixes such as '_mcon' or '_reamp' should not be included. Exclude additional identifiers that are specific to sample reps or recuts.
- **NOC**: #p
- **Conditioned**: condInitials (if applicable)
- **LR**: scenario for LR (refer to 3.5.3)
- **Case Notes** = a comment should be added here if a locus is ignored, indicating the reason: e.g. "D2S441 was ignored due to an unresolved allelic peak"

3.5.3 For LR scenarios, the naming format should start with the comparison sample's initials, followed by any conditioned samples' initials, and then the number and "U" for unknowns, followed by a "v" to separate the numerator hypothesis from the denominator hypothesis.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 6 OF 43

Examples	Resulting STRmix™ file name
Evidence File	
3-person, deconvolution, no conditioning, no comparisons	22-01234_567_1_1.1_trig_GS 3p 25-01234 1.1 3p
4-person deconvolution, no conditioning, comparing elimAB	22-01234_567_1_1.1_shirt_BL 4p AB3Uv4U 25-01234 1.1 4p AB3Uv4U
3-person deconvolution, conditioning vic CD, comparing elimAB	22-01234_567_1_1.1_shirt_BL 3p condCD-ABCD1UvCD2U 25-01234 1.1 3p condCD ABCD1UvCD2U
4-person deconvolution, conditioning elims CD, EF and GH, no comparisons	22-01234_567_1_1.1_trig_GS 4p condCDEFGH 25-01234 1.1 4p condCDEFGH
Suspect File	
1-person deconvolution, no conditioning, comparing suspTS	S22-05678 22-01234_567_2_1.1_slide_GS 1p TSv1U 25-05678 25-01234 1.1 1p TSv1U

3.6 Select , and confirm the settings against the following screenshot. They should always be the same for every STRmix™ analysis unless an exception is listed below. Any changes that are made will appear in bold on the run report.

- **MCMC Settings:** Burn-in Accepts (per chain) and Post Burn-in Accepts (per chain) must not be modified without documented approval from the technical leader (or his/her designee).
- **Mx Priors:** this will not be modified without documented approval from the technical leader (or his/her designee). See section 3.8 for more information regarding Mx priors.
- **Performance, Number of Threads:** it is okay to proceed if the Number of Threads does not match the screenshot below; this is specific to the computer being used.
- **Performance, Low Memory Mode:** This setting allows the computer to minimize the memory used for the run and can be turned on in order to use less computer memory (ex. if using your computer for other tasks while running STRmix™ in the background) or if a run fails to finish due to computer memory. Using this setting will increase the run time.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 7 OF 43

RUN SETTINGS

MCMC

Number of Chains 8	Burn-in Accepts (per chain) 10,000	Post Burn-in Accepts (per chain) 50,000
Random Walk SD 0.005	Post Burn-in Shortlist 9	☐ Extended Output

GELMAN-RUBIN

☐ Auto-Continue on GR
--

MX PRIORS

☐ Use Mx Priors
--

PERFORMANCE

Number of Threads 8	☐ Low Memory Mode
------------------------	--

SEED

☒ Random	724879
--	--------

- 3.7 If you are not adjusting Mx Priors settings, skip to 3.8; User informed Mx Priors is a function within STRmix™ that allows analysts to set approximate mixture proportion percentages for each contributor. If a proposed genotype does not fit the proportion percentage and variance set by the user, a penalty is applied to that iteration indicating an overall poor fit to the observed profile. Mx Priors settings are adjusted only with documented approval from the technical leader (or his/her designee). An analyst must show the Technical leader (or designee): 1) STRmix run data where a sample shows counterintuitive information in the output report and 2) the observed sample profile from the electropherogram. In addition, other troubleshooting techniques must be considered and tried prior to seeking approval for applying Mx Priors. See section 3.20 for set-up instructions for this function.
- 3.8 Select **Cancel** if all run settings are correct. If approved changes were made, Select **Apply**.
- 3.9 Set the number of contributors and ensure that contributor range remains unchecked. Select **Next**.
- 3.10 Check that the **Profiling Kit** selected is **Fusion_3500_2026**, **OCME_Fusion_3500** or **OCME_fusion_3130_2_7** depending on the data being analyzed. The STRmix™ output folder and file name (based on the combined Case Number and Sample ID you previously entered) are

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 8 OF 43

located at the top right. The kit labeled OCME_Fusion must not be used for deconvolutions in STRmix™ v2.7.

- 3.11 Click on **Kit Settings**. There are three tabs of settings to verify against the following screenshots in this window - **General**, **Stutters**, and **Loci**. If a locus needs to be ignored for a deconvolution, this is where you will be able to do that. See 0 for further instruction on ignoring a locus.

3.11.1 For the **Fusion_3500_2026 Profiling Kit**, ensure settings match the following:

KIT SETTINGS (Fusion_3500_2026)

GENERAL STUTTERS LOCI

VARIANCE

Allelic Variance 6.868, 3.785	Locus Amplification Variance 0.027	Minimum Variance Factor 0.5
Variance Minimization Parameter 1,000		

DROP-IN

Drop-in Cap 300	Drop-in Frequency 0.0087	Drop-in Distribution Parameters ☒ Uniform 22.31, 2.65
--------------------	-----------------------------	--

ADDITIONAL THRESHOLDS

Maximum Degradation 0.01	Degradation Start Point ☒ Use Smallest Peak	Saturation Threshold 30,000
-----------------------------	--	--------------------------------

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		
Document ID: 57029		
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 9 OF 43

KIT SETTINGS (Fusion_3500_2026)

GENERAL

STUTTERS

LOCI

BACK STUTTER

Maximum Stutter Ratio
☒ No Maximum 0.3

Variance
2.815, 1.174

FORWARD STUTTER

Maximum Stutter Ratio
☒ No Maximum 0.2

Variance
3.321, 3.202

HALF BACK (-2BP) STUTTER

Maximum Stutter Ratio
☒ No Maximum 0.1

Variance
1.645, 13.704

DOUBLE BACK (-8BP) STUTTER

Maximum Stutter Ratio
☒ No Maximum 0.1

Variance
3.051, 4.274

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 10 OF 43

3.11.1.1

KIT SETTINGS (Fusion_3500_2026)

GENERAL
STUTTERS
LOCI

LOCUS NAME	GENDER?	REPEAT LENGTH	IGNORE?	DETECTION THRESHOLD	BACK STUTTER	FORWARD STUTTER	HALF BACK (-2BP) STUTTER	DOUBLE BACK (-8BP) STUTTER
AMEL	☒							
D3S1358	☐	4	☐	85	☒	☒	☒	☒
D1S1656	☐	4	☐	85	☒	☒	☒	☒
D2S441	☐	4	☐	85	☒	☒	☒	☒
D10S1248	☐	4	☐	85	☒	☒	☒	☒
D13S317	☐	4	☐	85	☒	☒	☒	☒
Penta E	☐	5	☐	85	☒	☒	☒	☒
D16S539	☐	4	☐	120	☒	☒	☒	☒
D18S51	☐	4	☐	120	☒	☒	☒	☒
D2S1338	☐	4	☐	120	☒	☒	☒	☒
CSF1PO	☐	4	☐	120	☒	☒	☒	☒
Penta D	☐	5	☐	120	☒	☒	☐	☐
TH01	☐	4	☐	130	☒	☒	☒	☒
vWA	☐	4	☐	130	☒	☒	☒	☒
D21S11	☐	4	☐	130	☒	☒	☒	☒
D7S820	☐	4	☐	130	☒	☒	☒	☒
D5S818	☐	4	☐	130	☒	☒	☒	☐
TPOX	☐	4	☐	130	☒	☒	☒	☒
DYS391	☐	4	☒	130	☐	☐	☐	☐
D8S1179	☐	4	☐	160	☒	☒	☒	☒
D12S391	☐	4	☐	160	☒	☒	☒	☒
D19S433	☐	4	☐	160	☒	☒	☐	☒
FGA	☐	4	☐	160	☒	☒	☐	☒
D22S1045	☐	3	☐	160	☒	☒	☐	☒

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		
Document ID: 57029		
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 11 OF 43

3.11.2 For the **OCME_Fusion_3500 Profiling Kit**, ensure settings match the following:

KIT SETTINGS (OCME_Fusion_3500)

GENERAL

STUTTERS

LOCI

VARIANCE

Allelic Variance
7.427, 3.479

Locus Amplification Variance
0.03

Minimum Variance Factor
0.5

Variance Minimization Parameter
1,000

DROP-IN

Drop-in Cap
300

Drop-in Frequency
0.0087

Drop-in Distribution Parameters
☒ Uniform 22.31, 2.65

ADDITIONAL THRESHOLDS

Maximum Degradation
0.01

Degradation Start Point
☒ Use Smallest Peak

Saturation Threshold
30,000

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 12 OF 43

KIT SETTINGS (OCME_Fusion_3500)

GENERALSTUTTERSLOCITUTTERS

BACK STUTTER

Maximum Stutter Ratio

☐ No Maximum 0.3

Variance

1.799, 19.052

FORWARD STUTTER

Maximum Stutter Ratio

☐ No Maximum 0.2

Variance

1.999, 11.703

HALF BACK (-2BP) STUTTER

Maximum Stutter Ratio

☐ No Maximum 0.1

Variance

2.597, 6.632

DOUBLE BACK (-8BP) STUTTER

Maximum Stutter Ratio

☐ No Maximum 0.1

Variance

2.816, 6.792

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 13 OF 43

KIT SETTINGS (OCME_Fusion_3500)

GENERAL
STUTTERS
LOCI

LOCUS NAME	GENDER?	REPEAT LENGTH	IGNORE?	DETECTION THRESHOLD	BACK STUTTER	FORWARD STUTTER	HALF BACK (-2BP) STUTTER	DOUBLE BACK (-8BP) STUTTER
AMEL	☒							
D3S1358	☐	4	☐	85	☒	☒	☒	☐
D1S1656	☐	4	☐	85	☒	☒	☒	☒
D2S441	☐	4	☐	85	☒	☒	☒	☒
D10S1248	☐	4	☐	85	☒	☒	☒	☐
D13S317	☐	4	☐	85	☒	☒	☒	☐
Penta E	☐	5	☐	85	☒	☒	☐	☐
D16S539	☐	4	☐	120	☒	☒	☐	☐
D18S51	☐	4	☐	120	☒	☒	☒	☒
D2S1338	☐	4	☐	120	☒	☒	☒	☒
CSF1PO	☐	4	☐	120	☒	☐	☐	☐
Penta D	☐	5	☐	120	☒	☐	☐	☐
TH01	☐	4	☐	130	☒	☐	☒	☐
vwA	☐	4	☐	130	☒	☒	☐	☐
D21S11	☐	4	☐	130	☒	☒	☐	☒
D7S820	☐	4	☐	130	☒	☒	☒	☐
D5S818	☐	4	☐	130	☒	☒	☐	☐
TPOX	☐	4	☐	130	☒	☐	☐	☐
DYS391	☐	4	☒	130	☐	☐	☐	☐
D8S1179	☐	4	☐	160	☒	☒	☐	☒
D12S391	☐	4	☐	160	☒	☒	☐	☒
D19S433	☐	4	☐	160	☒	☒	☐	☐
FGA	☐	4	☐	160	☒	☒	☐	☒
D22S1045	☐	3	☐	160	☒	☒	☐	☐

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 14 OF 43

3.11.3 For the **OCME_Fusion_3130_2_7 Profiling Kit**, ensure settings match the following:

KIT SETTINGS (OCME_Fusion_3130)

GENERALSTUTTERSLOCI

VARIANCE

Allelic Variance
7.092, 0.933

Locus Amplification Variance
0.007

Minimum Variance Factor
0.5

Variance Minimization Parameter
1,000

DROP-IN

Drop-in Cap
100

Drop-in Frequency
0.0024

Drop-in Distribution Parameters
☒ Uniform

ADDITIONAL THRESHOLDS

Maximum Degradation
0.01

Degradation Start Point
☒ Use Smallest Peak

Saturation Threshold
8,000

KIT SETTINGS (OCME_Fusion_3130)

GENERALSTUTTERSLOCI

BACK STUTTER

Maximum Stutter Ratio
☐ No Maximum 0.3

Variance
2.108, 3.509

FORWARD STUTTER

Maximum Stutter Ratio
☐ No Maximum 0.1

Variance
2.312, 3.259

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		
Document ID: 57029		
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 15 OF 43

KIT SETTINGS (OCME_Fusion_3130)

GENERAL
STUTTERS
LOCI

LOCUS NAME	GENDER?	REPEAT LENGTH	IGNORE?	DETECTION THRESHOLD	BACK STUTTER	FORWARD STUTTER
AMEL	☒					
D3S1358	☐	4	☐	50	☒	☒
D1S1656	☐	4	☐	50	☒	☒
D2S441	☐	4	☐	50	☒	☒
D10S1248	☐	4	☐	50	☒	☒
D13S317	☐	4	☐	50	☒	☒
Penta E	☐	5	☐	50	☒	☒
D16S539	☐	4	☐	50	☒	☒
D18S51	☐	4	☐	50	☒	☒
D2S1338	☐	4	☐	50	☒	☒
CSF1PO	☐	4	☐	50	☒	☒
Penta D	☐	5	☐	50	☒	☒
TH01	☐	4	☐	50	☒	☒
vWA	☐	4	☐	50	☒	☒
D21S11	☐	4	☐	50	☒	☒
D7S820	☐	4	☐	50	☒	☒
D5S818	☐	4	☐	50	☒	☒
TPOX	☐	4	☐	50	☒	☒
DYS391	☐	4	☒	50	☒	☒
D8S1179	☐	4	☐	50	☒	☒
D12S391	☐	4	☐	50	☒	☒
D19S433	☐	4	☐	50	☒	☒
FGA	☐	4	☐	50	☒	☒
D22S1045	☐	3	☐	50	☒	☒

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 16 OF 43

- 3.12 Once all settings have been verified, click **Cancel** to return to the previous screen. For any **locations that need to be ignored**, click the checkbox in the **Ignore** column before clicking **Apply** and make sure a note is made in the **Case Notes** section of the interpretation.
- 3.13 Add your **Evidence Profile Data** by either dragging and dropping the .txt file into the box or clicking on the **+** to find your .txt file. The following window will pop up and you can select your evidence profile(s). If multiple replicates are in the same run, you can select multiple. Click **Confirm** to return to the file import screen.
- 3.13.1 Repeat for any replicate(s) from other .txt files.

ADD EVIDENCE PROFILE DATA (OCME_Fusion_3130)

YY-XXXX - OCME ID #p condABCD (if applicable) HPvHD (if applicable)

Input File Carmody010222 44-47U A_SM.txt Select

SAMPLE NAME

- ☐ 02-Amp_Pos_12282021.141254_B01_04.hid
- ☐ 03-Amp_Neg_12282021.141254_C01_07.hid
- ☐ 04-E_Neg1_12132021.151658_D01_10.hid
- ☐ 05-FB21-07267_116_2_1.1_strap1_GS_E01_13.hid
- ☐ 06-FB21-07267_116_2_1.1_strap1_GS_rep_F01_16.hid
- ☒ 07-FB21-07267_116_4_1.1_trigg2_GS_G01_19.hid
- ☒ 08-FB21-07267_116_4_1.1_trigg2_GS_rep_H01_22.hid
- ☐ 09-FB21-07267_116_5_1.1_strap2_GS_A02_02.hid
- ☐ 10-FB21-07267_116_5_1.1_strap2_GS_rep_B02_05.hid
- ☐ 11-FB21-07267_116_6_1.1_slide_GS_C02_08.hid
- ☐ 12-FB21-07258_778_2_1.1_straps_GS_D02_11.hid
- ☐ 13-FB21-07258_778_2_1.1_straps_GS_rep_E02_14.hid
- ☐ 14-FB21-07258_778_5_1.1_straps_GS_F02_17.hid
- ☐ 15-FB21-07258_778_5_1.1_straps_GS_rep_G02_20.hid

☒ ☐

Cancel Confirm

- 3.14 To add a conditioned contributor, add your **Reference Profile Data** by either dragging and dropping the .txt file into the box or clicking on the **+** to find your .txt file. If you are performing a deconvolution without a conditioned contributor, proceed to step 3.15.
- 3.14.1 Select your reference profile(s) for conditioning using the checkboxes and click **Confirm** (similar to the screen for selecting your evidence profile).

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 17 OF 43

- 3.14.2 Once the profile data is added, the box for HP will be checked on the right side. To condition, you must also check the box for HD. Conditioned contributors are considered true donors in Hp and Hd.

REFERENCE PROFILE DATA

+

13-FBPT22-00001_22524_1_1_victim_E02_14.hid_REF.csv

CONTRIBUTOR TO:

HP

☒

HD

☒

- 3.15 To perform an **LR in conjunction with the deconvolution**, add your **Reference Profile Data** by either dragging and dropping the .txt file into the box or clicking on the + to find your .txt file. If you are performing a deconvolution without an LR, proceed to step 3.16.

- 3.15.1 Select your reference profile for LR's using the checkboxes and click Confirm (similar to the screen for selecting your evidence profile). The reference will be added within the numerator of the LR (i.e. assigned to HP only).

- 3.15.2 Once the reference profile is added to your HP, you can then select LR Settings, and check them against the following screenshot:

LR SETTINGS

POPULATIONS

NAME	PROPORTION	FST	ALLELE FREQUENCY FILE
nist_afam_2_7	0.25	0.03b(1.0, 1.0)	NIST Fusion AfAm_Amended20...
nist_asian_2_7	0.25	0.03b(1.0, 1.0)	NIST Fusion Asian_Amended20...
nist_cauc_2_7	0.25	0.03b(1.0, 1.0)	NIST Fusion Cauc_Amended201...
nist_hisp_2_7	0.25	0.03b(1.0, 1.0)	NIST Fusion Hisp_Amended201...

+

-

📁

SUB-SOURCE LR

☒ Assign Sub-Source LR

SAMPLING VARIATION

☒ Calculate HPD

Number of HPD Iterations
1000

☒ MCMC Uncertainty

Probability Interval Quantile
99

☒ Allele Frequency Uncertainty

Probability Interval Sides
1


FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 18 OF 43

3.15.3 Select  to return to the file import screen.

3.16 Select  (or  for Batch Mode, Section 6). The Progress window will open showing the Pre-Burnin, Burnin progress, and MCMC Progress (or for Batch Mode, you will return to the batch screen to continue to add to the batch).

3.16.1 If a flag fires about missing an expected stutter (see below), review the input data and confirm that an allele label was not inadvertently deleted during analysis and that the input .txt file is correct before proceeding.

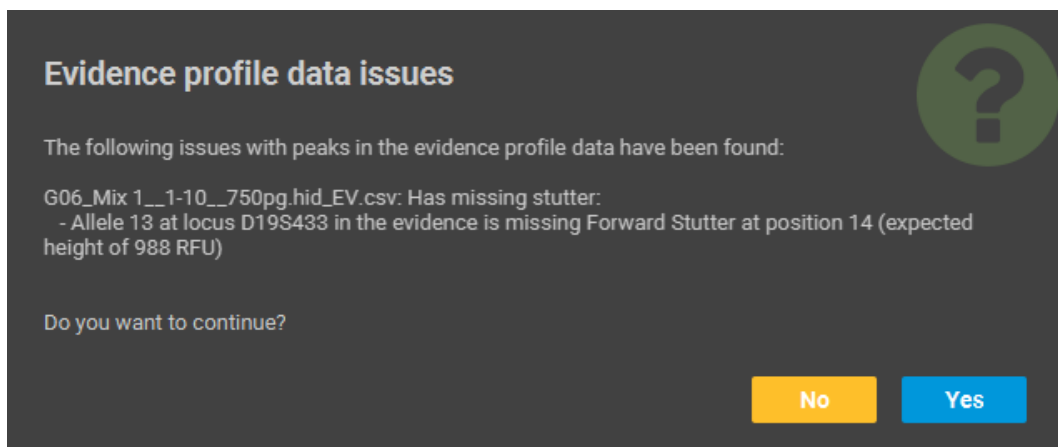


3.16.1.1 If no editing changes need to be made, you can click **Yes** to continue. If changes need to be made, select **No**, correct the labeling in the electropherogram and see Section 2.1.5.6 for editing the input file. Return to the beginning of setting up the run.

3.16.2 If a flag fires about peaks below analytical threshold (see below). Click  and do not continue. Review the input data and the GeneMarker run settings used for analysis.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 19 OF 43

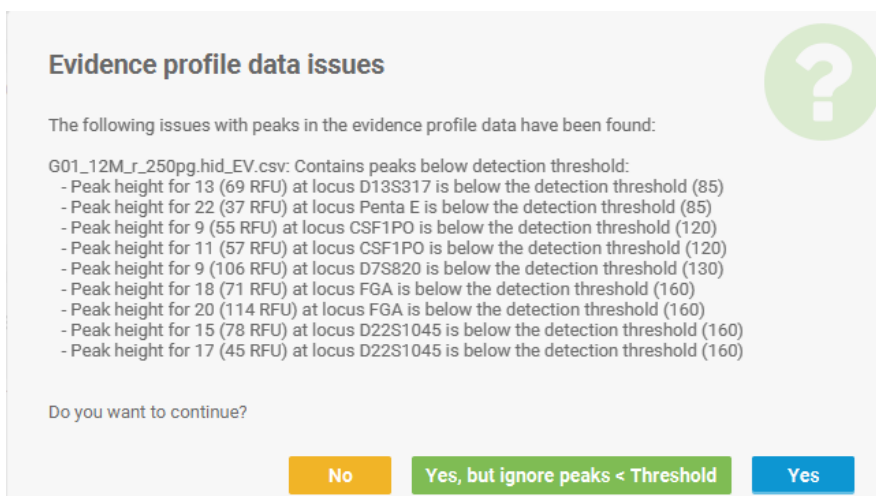

3.17 Once complete, you will be at the Results screen.

3.17.1 You can click on the single file folder icon  on the right side to be directed to that run's **Results** folder. Select to return to the Main Menu.

3.17.2 Alternatively, you can manually navigate to your STRmix™ Run Folder(s) within the **OCME_STRmix_Fileshare** folder (\\csc\ocme\OCME_STRmix_Fileshare).

3.18 **COPY** your run folder(s) into the previously created FB sample folder within the **STRmix Data** folder.

3.19 Once you have copied the folder, **CONFIRM that all files for that run have transferred over correctly** to the **STRmix Data** folder. After confirmation, the copy of the STRmix™ Run folder located in the **OCME_STRmix_Fileshare** folder should be **deleted**.

3.20 When Applying User Informed Mx Priors: ****DOCUMENTED APPROVAL BY TECHNICAL LEADER NEEDED PRIOR TO APPLICATION****


3.20.1 Use the stutter filtered electropherograms to determine an approximate mixture proportion percentage for each contributor. Start within the first few loci to the left in each dye channel.

3.20.2 Set-up an interpretation run as per section 3.1-3.6 of this manual. Be sure your “Number of Contributors” setting is correct, and the “Contributor Range” option remains unchecked.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 20 OF 43

- 3.20.3 Select **Run Settings** to open the RUN SETTINGS window and check the box marked “Use Mx Priors”. A CONFIGURE MX PRIORS window will appear like the one seen below:

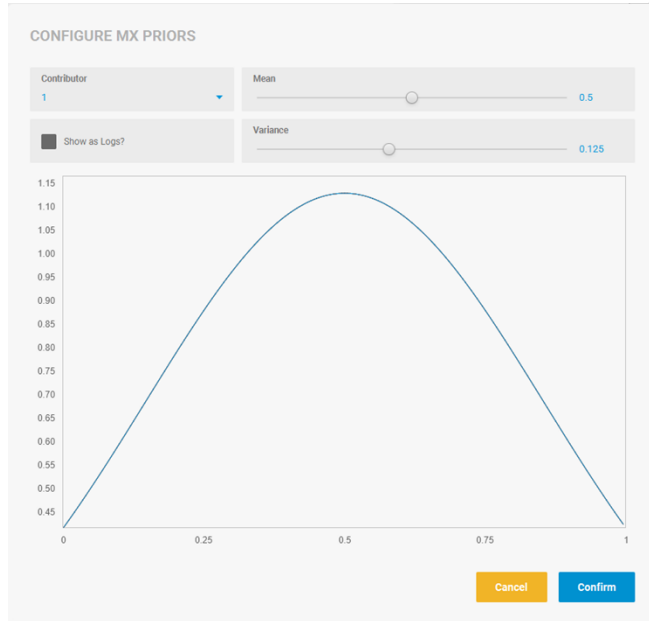


- 3.20.4 Ensure contributor 1 is chosen in the drop-down menu. Using the top slider bar marked “Mean” to set the contributor proportion to the previously determined value.
- 3.20.5 Look at the range of the mixture proportion values used to determine the mean contributor as a guide for setting the bottom slider bar marked “Variance”. Suggestions for variance to be set based on previous testing range between 0.125 and 6.1035E-5.
- 3.20.6 Repeat the steps in sections 3.20.4 and 3.20.5, changing the mean setting for the additional contributor(s) in the sample based on the determined values as needed. The variance setting should be set to the same value for all contributors.
- 3.20.7 Tap **Confirm** in the CONFIGURE MX PRIORS window, then click **Apply** in the RUN SETTINGS window.
- 3.20.8 Click **Next**.
- 3.20.9 Go back to section 3.10 and follow the steps needed to complete the run.

4 Likelihood Ratio Calculations with STRmix™ (Investigation: LR from Previous)

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 21 OF 43

- 4.1 Note: This section is specific for running LR's using a previously run deconvolution. Samples must undergo deconvolution prior to (or in conjunction with) running an LR for a comparison sample. Refer to [Section 3 for setting up deconvolutions in conjunction with the LR](#) and the [Interpretation of PowerPlex® Fusion data run on 3500xL](#) manual for further guidance on running deconvolutions and LR's.
- 4.2 From the **Main Menu**, select **Investigation**. Within the **Investigation** screen, select **LR from Previous**. You will then select the deconvolution file for which you would like to calculate an LR for a comparison sample:

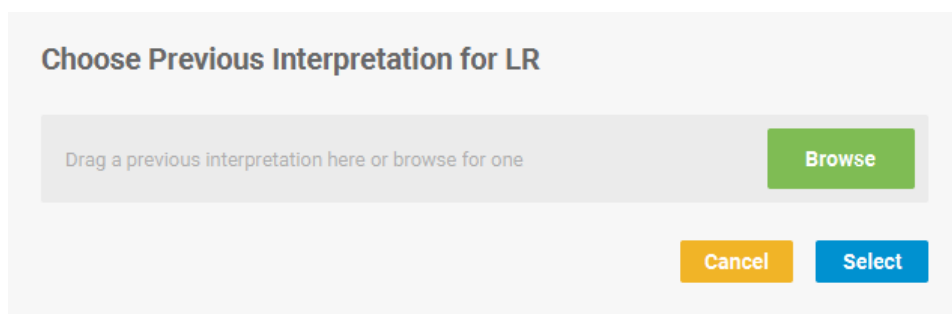

- 4.3 Drag and drop the entire decon run folder or select Browse to locate **the deconvolution file**, labeled as config.xml or results.xml (or settings.ini for older versions of STRmix™) within your run folder in the STRmix Data folder in the M drive. Click **Select** to navigate to the **LR from Previous** screen to name the run.
- 4.4 When naming STRmix™ Likelihood Ratio runs, the Case Number and Sample ID will auto-populate from the deconvolution file. These should be updated to the appropriate naming convention for an LR run before proceeding.
- 4.5 The following naming convention should be used for likelihood ratios:

4.5.1 Evidence Files

- **Case Number** = YY-XXXXX (do not include "FB" or "F")
- **Sample ID** = item and sample identifiers consistent with evidence electropherogram. Suffixes such as '_mcon' or '_reamp' should not be included. Exclude additional identifiers that are specific to sample reps or recuts.
- **NOC**: #p
- **Conditioned**: condInitials (if applicable)
- **LR**: scenario for LR (if applicable, refer to 4.5.3)
- **Case notes** = a comment should be added here if a locus is ignored, indicating the reason: e.g. "D2S441 was ignored due to an unresolved allelic peak"

4.5.2 Suspect Files

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 22 OF 43

- **Case Number** = SYX-XXXXX (do not include “FB” or “F”)
- **Sample ID** = item and sample identifiers consistent with evidence electropherogram. (include evidence file FB# without the “FB”) Suffixes such as ‘_mcon’ or ‘_reamp’ should not be included. Exclude additional identifiers that are specific to sample reps or recuts.
- **NOC**: #p
- **Conditioned**: condInitials (if applicable)
- **LR**: scenario for LR (refer to 4.5.3)
- **Case Notes** = a comment should be added here if a locus is ignored, indicating the reason: e.g. “D2S441 was ignored due to an unresolved allelic peak”

4.5.3 For LR scenarios, the naming format should start with the comparison sample’s initials, followed by any conditioned samples’ initials, and then the number and “U” for unknowns, followed by a “v” to separate the numerator hypothesis from the denominator hypothesis.

Examples	Resulting STRmix file name
Evidence File	
comparing elimAB to 4p decon, no cond	22-01234_567_1_1.1_shirt_BL AB3Uv4U 25-01234_1.1 AB3Uv4U
comparing elimAB to 2p decon, cond vicCD	22-01234_567_1_1.1_shirt_BL ABCDvCD1U 25-01234_1.1 ABCDvCD1U
comparing elimAB to 3p decon, cond vicCD and elimEF	22-01234_567_1_1.1_shirt_BL ABCDEFvCDEF1U 25-01234_1.1 ABCDEFvCDEF1U
Suspect File	
comparing suspTS to 1p decon, no cond	S22-05678_22-01234_567_2_1.1_slide_GS TSv1U 25-05678_25-01234_1.1 TSv1U
comparing suspTS to 2p decon, no cond	S22-05678_22-01234_567_2_1.1_slide_GS TS1Uv2U 25-05678_25-01234_1.1 TS1Uv2U
comparing suspTS to 3p decon, cond vicCD	S22-05678_22-01234_567_2_1.1_slide_GS TSCD1UvCD2U 25-05678_25-01234_1.1 TSCD1UvCD2U
comparing suspTS to 4p decon, cond vicCD and elimEF	S22-05678_22-01234_567_2_1.1_slide_GS TSCDEF1UvCDEF2U 25-05678_25-01234_1.1 TSCDEF1UvCDEF2U

4.6 The run settings will be pulled from the deconvolution and should not be changed. To check, click on **Run Settings** and see 3.6 for Run Settings screenshot. Click **Cancel** if all settings are correct and click Next to proceed.

- **Performance, Number of Threads**: it is okay to proceed if the Number of Threads is different; this is specific to the computer being used.
- **Performance, Low Memory Mode**: This setting allows the computer to minimize the memory used for the run and can be turned on in order to use less computer memory (ex. if using your computer for other tasks while running STRmix™ in the background) or if a run fails to finish due to computer memory. Using this setting will increase the run time.
- if a run fails to finish due to computer memory.

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 23 OF 43

- **Performance, Seed:** This number may be different; there is a Seed for each run, deconvolution or LR.
- 4.7 The kit settings will also be pulled from the deconvolution. The kit will be pulled in as appropriate based on the STRmix deconvolution file:
- 4.7.1 **Fusion_3500_2026** or OCME_Fusion_3500 for 3500xL evidence data.
- 4.7.2 OCME_Fusion_3130_2_7 for 3130xL evidence data deconvoluted using STRmix v2.7.
- 4.7.3 **OCME_Fusion** for 3130xL evidence data deconvoluted using STRmix v2.4.
- 4.8 If a locus needs to be ignored, click on the Kit settings.
- 4.8.1 Rarely, a locus may be ignored at this step. For example, in the case of a tri-allelic pattern that matches your reference sample, which was not recognized at the deconvolution stage. For any **locations that need to be ignored**, click the checkbox in the **Ignore** column before clicking .
- 4.9 Your Evidence Profile Data and Reference Profile Data for any conditioned contributors (HP & HD selected) will be populated from the deconvolution.
- 4.10 Add the **profile(s) for comparison** to the **Reference Profile Data** by either dragging and dropping the .txt file into the box or clicking on the  to find your .txt file(s).
- 4.10.1 Select your reference profile(s) for the LR using the checkboxes and click . The reference will be added within the numerator of the LR (i.e. assigned to HP only).
- 4.10.2 Once the reference profile(s) is added to your HP, you can then select , and check them against the following screenshot:

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 24 OF 43

LR SETTINGS

POPULATIONS

NAME	PROPORTION	FST	ALLELE FREQUENCY FILE
nist_afam_2_7	0.25	0.03b(1.0, 1.0)	NIST Fusion AfAm_Amended20...
nist_asian_2_7	0.25	0.03b(1.0, 1.0)	NIST Fusion Asian_Amended20...
nist_cauc_2_7	0.25	0.03b(1.0, 1.0)	NIST Fusion Cauc_Amended201...
nist_hisp_2_7	0.25	0.03b(1.0, 1.0)	NIST Fusion Hisp_Amended201...

+ - 

SUB-SOURCE LR

☒ Assign Sub-Source LR

SAMPLING VARIATION

☒ Calculate HPD ☒ MCMC Uncertainty ☒ Allele Frequency Uncertainty

Number of HPD Iterations: 1000 Probability Interval Quantile: 99 Probability Interval Sides: 1

- 4.11 Select **Start** (or **Queue** for Batch Mode, Section 6). The Progress window will open (or for Batch Mode, you will return to the batch screen to continue to add to the batch).
- 4.12 Once complete, you will be at the Results screen.
 - 4.12.1 You can click on the single file folder icon  on the right side to be directed to that run's **Results** folder. Select **Finish** to return to the Investigation Menu.
 - 4.12.2 Alternatively, you can manually navigate to your STRmix Run Folder within the **OCME_STRmix_Fileshare** folder (\\csc\ocme\OCME_STRmix_Fileshare).
- 4.13 **COPY** your run folder into the previously created FB sample folder within the **STRmix Data** folder.
- 4.14 Once you have copied the folder, **CONFIRM that all files for that run have transferred over correctly** to the **STRmix Data** folder. After confirmation, the copy of the STRmix Run folder located in the **OCME_STRmix_Fileshare** folder should be **deleted**.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 25 OF 43

5 LRs for Convicted Offender Match Testimony


- 5.1 For testimony in relation to a convicted offender match, where a statistic is needed, an LR must be calculated through STRmix™.
 - 5.1.1 Obtain the Convicted Offender profile through the CODIS software:
 - 5.1.1.1 Go to the Specimen Manager window → Click the filter button at the top of the Specimen ID column → Click “custom” → Enter your specimen ID from your evidence sample in the top line in the dialog box → Hit “search”
 - 5.1.1.2 Once your specimen ID pops up, right click on the line and click “view matches” → Right-click and choose “print reports”, choose “match details short report”
 - 5.1.2 Add this to your case file and use the convicted offender profile from this report in order to create a comparison sample profile for STRmix™ LR calculation using the form: [Reference Profile for STRmix Input Creation](#).
 - 5.1.3 Once the STRmix™ LR report is generated, it will need to undergo technical review (documented using case contacts) and recertification prior to testimony.

6 Running STRmix™ using Batch Mode

- 6.1 Several STRmix™ runs can be set up and queued to run sequentially. To set up a queued analysis for multiple runs, select **Batch Mode** from the STRmix™ main window. Batch mode should only be used for samples from one individual case within a single batch.
- 6.2 Before setting up a batch, navigate to the OCME_STRmix_Fileshare folder ([\\csc\ocme\OCME_STRmix_Fileshare](#)). Create a new folder within the STRmix Fileshare folder with the name of “Batch Mode [date] [your initials].”
- 6.3 Select  to change the Batch Directory. Navigate to and select the created folder inside the **OCME_STRmix_Fileshare**. This new folder will now appear at the top of the screen:

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 26 OF 43

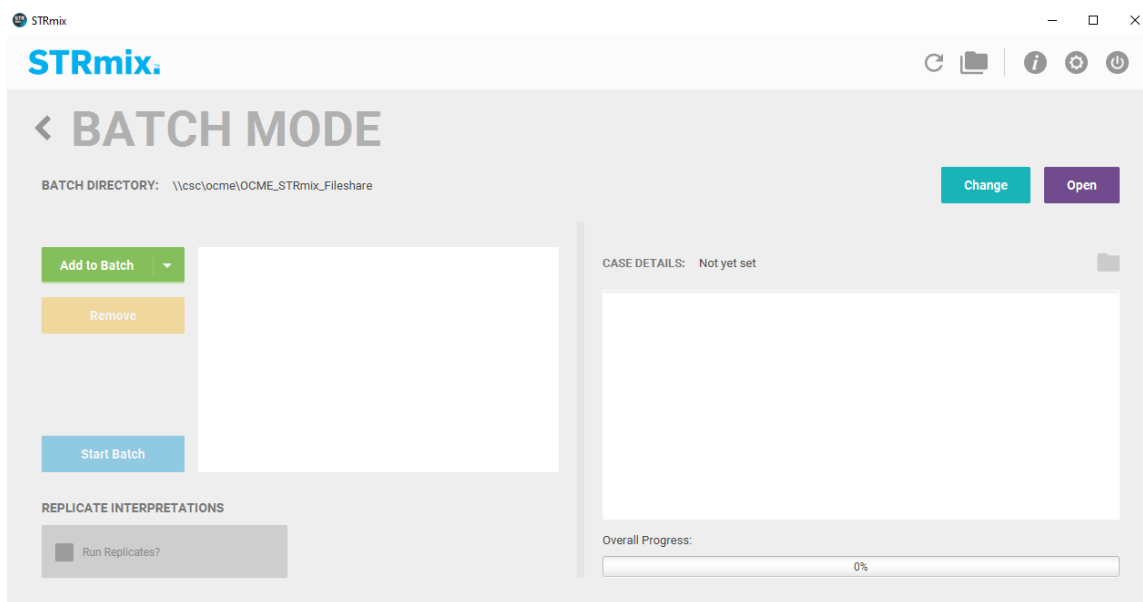

- 6.4 Use the down arrow next to **Add to Batch**  and **Add Interpretation** or **Add LR from Previous** to set up each individual run. If you do not click the arrow and click **Add to Batch**, it will automatically open whichever run type was used last.
 - 6.4.1 These batch details will be saved even if you close STRmix™ until you set up a new run; the batch can be stopped and restarted at a later date.
- 6.5 Complete the analysis set up for the first sample following the corresponding setup instructions: [Section 3 for Deconvolutions \(Investigation\)](#) or [Section 4 for LRs \(Investigation: LR from Previous\)](#).
- 6.6 After setting up your run, click  to return to the Batch Mode setup window.
- 6.7 Repeat steps 6.4-6.6 to add additional runs.
 - 6.7.1 If you need to edit a run, you will need to remove it and redo the setup. To remove a sample from the batch mode, highlight the case/sample in the **Calculations in Batch** part of the window and select **Remove**.
- 6.8 Select  to start the batch run.
- 6.9 After completion of analyses, use the left arrow at the top left next to Batch Mode to return to the STRmix™ Main Menu. You can click  to be directed to the Batch Mode folder you created.
 - 6.9.1 Alternatively, if you have additional runs to perform within the same batch, you can add additional runs and continue the batch.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 27 OF 43

- 6.9.2 The **Batch Log** will show if any of the individual runs failed and why. If an individual run failed, the run folder will not contain a Results folder, and a .txt file called BATCH CALCULATION FAILED will be made instead.
- 6.10 Results folders from Batch Mode will be saved in the folder you created within the **OCME_STRmix_Fileshare** folder ([\\csc\ocme\OCME_STRmix_Fileshare](#)).
- 6.11 **COPY your run folders** into the previously created FB sample folder(s) within the **STRmix Data** folder.
- 6.12 Once you have copied the folders, **CONFIRM that all files have transferred over correctly** to the **STRmix Data** folders. After confirmation, the copy of the STRmix Run folder and the Batch folder with your initials and date located in the **OCME_STRmix_Fileshare** folder should be **deleted**.

7 Evaluation of the STRmix™ Analysis Setup

- 7.1 Verify the evidence and reference input sections of the STRmix™ printout against the associated electropherograms. Ensure that:
- 7.1.1 All appropriate edits were made; no artifact peaks were left labeled and no allelic peak labels were removed.
 - 7.1.1.1 If an Evidence profile issue flag was generated at the beginning of your run and checked for any issues against the electropherogram prior to running, this will also show up in the **EVIDENCE PEAK ISSUES** section at the end of the report. The information listed here should be aligned with the issue flag.
 - 7.1.2 Correct input file(s) have been selected.
 - 7.1.3 The correct file was imported into an LR from previous analysis, if applicable.
 - 7.1.4 All suitable replicates have been utilized.
 - 7.1.4.1 Check the **Inter replicate efficiency** that is generated by STRmix™. If there are drastic inconsistencies between the two efficiencies that are consistent with the amount of data present in the replicates, the STRmix analysis may be marked as ‘not reported’ and a new analysis may be performed with the amplification containing the most information. A third amplification may be warranted.
- 7.2 The number of contributors that best describes the data was chosen, as applicable.
- 7.3 The correct assumptions (i.e. conditioning) have been made, if applicable.
- 7.4 The appropriate proposition has been selected (i.e. LR calculation), if applicable.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 28 OF 43

- 7.5 The **SEED** value listed at the beginning of the report, with the CASE NUMBER, SAMPLE NAME, and COMMENTS, is the starting number used within the random number generator.
- 7.6 Check the **SETTINGS** at the end of the report to verify that the STRmix™ run was set up properly; note that any edited settings values are bolded by the program.
- 7.6.1 The **CASE SETTINGS** will be specific to your case and should be used to check for correct setup with the run name.
- 7.6.2 The **MCMC SETTINGS** are only for runs that include a deconvolution.
- 7.6.2.1 Check the remaining settings against the following screenshot keeping in mind the numbers of accepts may be different (and bolded) if this setting has been approved for use.

MCMC SETTINGS

Number of contributors	4
Use Mx priors	N
Number of chains	8
Burn-in accepts per chain	10,000
Post burn-in accepts per chain	50,000
Random walk SD	0.005
Post burn-in shortlist	9.0
Auto-continue on Gelman-Rubin	N

- 7.6.3 The **KIT SETTINGS** will be included in deconvolutions and LR_s. If any loci were ignored, they will be listed here.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 29 OF 43

Fusion 3500 2026 Kit

KIT SETTINGS

Ignored loci	DYS391
Detection thresholds	D3S1358 85
	D1S1656 85
	D2S441 85
	D10S1248 85
	D13S317 85
	Penta E 85
	D16S539 120
	D18S51 120
	D2S1338 120
	CSF1PO 120
	Penta D 120
	TH01 130
	vWA 130
	D21S11 130
	D7S820 130
	D5S818 130
	TPOX 130
	D8S1179 160
	D12S391 160
	D19S433 160
	FGA 160
	D22S1045 160
Saturation	30,000
Degradation starts at	-1.0
Degradation max	0.01
Drop-in cap	300
Drop-in rate parameter	0.0087
Drop-in parameters (α , β)	22.31, 2.65
Min variance factor	0.5
Variance minimization parameter	1,000
LSAE variance parameter ($1/\lambda$)	0.027
Allelic variance parameters (α , β)	6.868, 3.785

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 30 OF 43

Back Stutter

Maximum stutter ratio

0.3

Stutter variance parameters (α , β)

2.815, 1.174

Forward Stutter

Maximum stutter ratio

0.2

Stutter variance parameters (α , β)

3.321, 3.202

Half Back (-2bp) Stutter

Maximum stutter ratio

0.1

Stutter variance parameters (α , β)

1.645, 13.704

Double Back (-8bp) Stutter

Maximum stutter ratio

0.1

Stutter variance parameters (α , β)

3.051, 4.274

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 31 OF 43

OCME_Fusion_3500 Kit

KIT SETTINGS

Ignored loci	DYS391	
Detection thresholds	D3S1358	85
	D1S1656	85
	D2S441	85
	D10S1248	85
	D13S317	85
	Penta E	85
	D16S539	120
	D18S51	120
	D2S1338	120
	CSF1PO	120
	Penta D	120
	TH01	130
	vWA	130
	D21S11	130
	D7S820	130
	D5S818	130
	TPOX	130
	D8S1179	160
	D12S391	160
	D19S433	160
	FGA	160
	D22S1045	160
Saturation	30,000	
Degradation starts at	-1.0	
Degradation max	0.01	
Drop-in cap	300	
Drop-in rate parameter	0.0087	
Drop-in parameters (α , β)	22.31, 2.65	
Min variance factor	0.5	
Variance minimization parameter	1,000	
LSAE variance parameter ($1/\lambda$)	0.03	
Allelic variance parameters (α , β)	7.427, 3.479	

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 32 OF 43

Back Stutter

Maximum stutter ratio	0.3
Stutter variance parameters (α , β)	1.799, 19.052

Forward Stutter

Maximum stutter ratio	0.2
Stutter variance parameters (α , β)	1.999, 11.703

Half Back (-2bp) Stutter

Maximum stutter ratio	0.1
Stutter variance parameters (α , β)	2.597, 6.632

Double Back (-8bp) Stutter

Maximum stutter ratio	0.1
Stutter variance parameters (α , β)	2.816, 6.792

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 33 OF 43

OCME_Fusion_3130_2_7 Kit

KIT SETTINGS

Ignored loci	DYS391
Detection thresholds	D3S1358 50
	D1S1656 50
	D2S441 50
	D10S1248 50
	D13S317 50
	Penta E 50
	D16S539 50
	D18S51 50
	D2S1338 50
	CSF1PO 50
	Penta D 50
	TH01 50
	vWA 50
	D21S11 50
	D7S820 50
	D5S818 50
	TPOX 50
	D8S1179 50
	D12S391 50
	D19S433 50
	FGA 50
	D22S1045 50
Saturation	8,000
Degradation starts at	-1.0
Degradation max	0.01
Drop-in cap	100
Drop-in rate parameter	0.0024
Drop-in parameters (α , β)	0.0, 0.0
Min variance factor	0.5
Variance minimization parameter	1,000
LSAE variance parameter ($1/\lambda$)	0.007
Allelic variance parameters (α , β)	7.092, 0.933
Back Stutter	
Maximum stutter ratio	0.3
Stutter variance parameters (α , β)	2.108, 3.509
Forward Stutter	
Maximum stutter ratio	0.1
Stutter variance parameters (α , β)	2.312, 3.259

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 34 OF 43

OCME_Fusion Kit (for likelihood ratios on 3130xL STRmix™ v2.4 data only)

KIT SETTINGS

Ignored loci	DYS391
Detection thresholds	D3S1358 50
	D1S1656 50
	D2S441 50
	D10S1248 50
	D13S317 50
	Penta E 50
	D16S539 50
	D18S51 50
	D2S1338 50
	CSF1PO 50
	Penta D 50
	TH01 50
	vWA 50
	D21S11 50
	D7S820 50
	D5S818 50
	TPOX 50
	D8S1179 50
	D12S391 50
	D19S433 50
	FGA 50
	D22S1045 50
Saturation	8,000
Degradation starts at	-1.0
Degradation max	0.01
Drop-in cap	100
Drop-in rate parameter	0.0024
Drop-in parameters (α , β)	0.0, 0.0
Min variance factor	0.5
Variance minimization parameter	1,000
LSAE variance parameter ($1/\lambda$)	0.0065
Allelic variance parameters (α , β)	9.1374, 0.7472
<i>Back Stutter</i>	
Maximum stutter ratio	0.3
Stutter variance parameters (α , β)	1.5007, 12.9748
<i>Forward Stutter</i>	
Maximum stutter ratio	0.1
Stutter variance parameters (α , β)	9.1374, 0.7472

- 7.6.4 The **PROFILE SETTINGS** can be used to check input file names and that the run was set up properly. If a stutter flag was generated as seen in step 3.16, “Ignore peaks below detection threshold” may be present as seen below; however this does not indicate that there are peaks below threshold within the run.

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 35 OF 43

PROFILE SETTINGS

Number of evidence profiles	1
Evidence profile filenames	63-FB21-07367_658_1.4.1_1.1_DS1_ST_G08_20.hid_EV.csv
Number of Hp knowns	0
Number of Hd knowns	0
Ignore peaks below detection threshold	N

7.6.5 The **LR SETTINGS** will only be present when an LR was run (whether separately or along with the deconvolution).

LR SETTINGS

Number of populations	4
Assign sub-source LR	Y
Calculate HPD	Y
HPD iterations	1,000
Use MCMC uncertainty	Y
Use Allele Frequency uncertainty	Y
HPD quantile	99%
HPD sides	1
NIST_AFAM_2.7	
Proportion	0.25
FST	0.03b(1.0, 1.0)
Allele frequency file	NIST Fusion AfAm_Amended2017.csv
NIST_ASIAN_2.7	
Proportion	0.25
FST	0.03b(1.0, 1.0)
Allele frequency file	NIST Fusion Asian_Amended2017.csv
NIST_CAUC_2.7	
Proportion	0.25
FST	0.03b(1.0, 1.0)
Allele frequency file	NIST Fusion Cauc_Amended2017.csv
NIST_HISP_2.7	
Proportion	0.25
FST	0.03b(1.0, 1.0)
Allele frequency file	NIST Fusion Hisp_Amended2017.csv

7.6.6 The **PERFORMANCE SETTINGS** shows the number of threads (specific to the computer used to run) and whether Low Memory Mode was used.

7.6.7 The **ADDITIONAL KIT DETAILS** lists the stutter files used in a deconvolution analysis.

Controlled versions of Department of Forensic Biology Manuals only exist in the Forensic Biology Qualtrax software. All printed versions are non-controlled copies.

© NYC OFFICE OF CHIEF MEDICAL EXAMINER

Qualtrax template 040621

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 36 OF 43

Fusion 3500 2026 Kit

ADDITIONAL KIT DETAILS

Size Regression File	Fusion_SizeRegression.csv
Back Stutter	
Position Relative to Parent	-1, 0
Inversely Proportional To	Expected Height of Stutter Peak
Stutter Regression File	OCME_BSFile_101325.txt
Stutter Exceptions File	OCME_BSEExceptionsFile_101325.csv
Forward Stutter	
Position Relative to Parent	1, 0
Inversely Proportional To	Expected Height of Stutter Peak
Stutter Regression File	OCME_FSFile_101325.txt
Stutter Exceptions File	OCME_FSEExceptionsFile_101325_None.csv
Half Back (-2bp) Stutter	
Position Relative to Parent	0, -2
Inversely Proportional To	Expected Height of Stutter Peak
Stutter Regression File	OCME_HBSFile_101325.txt
Stutter Exceptions File	OCME_HBSEExceptionsFile_101325_None.csv
Double Back (-8bp) Stutter	
Position Relative to Parent	-2, 0
Inversely Proportional To	Expected Height of Stutter Peak
Stutter Regression File	OCME_DBSFile_101325.txt
Stutter Exceptions File	OCME_DBSEExceptionsFile_101325_None.csv

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 37 OF 43

OCME Fusion 3500 Kit ADDITIONAL KIT DETAILS

Size Regression File	fusion_sizeregression.csv
Back Stutter	
Position Relative to Parent	-1, 0
Inversely Proportional To	Observed Height of Parent Allele
Stutter Regression File	OCME_BackStutterFile3500_022321.txt
Stutter Exceptions File	OCME_BSEExceptionsFile_020321.csv
Forward Stutter	
Position Relative to Parent	1, 0
Inversely Proportional To	Expected Height of Stutter Peak
Stutter Regression File	OCME_ForwardStutterFile3500_022321.txt
Stutter Exceptions File	OCME_FSEExceptions File_092020.csv
Half Back (-2bp) Stutter	
Position Relative to Parent	0, -2
Inversely Proportional To	Expected Height of Stutter Peak
Stutter Regression File	OCME_HalfBackStutterFile3500_022321.txt
Stutter Exceptions File	OCME_HBSEExceptionsFile_092020.csv
Double Back (-8bp) Stutter	
Position Relative to Parent	-2, 0
Inversely Proportional To	Expected Height of Stutter Peak
Stutter Regression File	OCME_DoubleBackStutterFile3500_022321.txt
Stutter Exceptions File	OCME_DBSEExceptionsFile_092020.csv

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 38 OF 43

OCME_Fusion_3130_2_7 Kit

ADDITIONAL KIT DETAILS

Size Regression File	Fusion_SizeRegression.csv
<i>Back Stutter</i>	
Position Relative to Parent	-1, 0
Inversely Proportional To	Observed Height of Parent Allele
Stutter Regression File	OCME_Fusion_Stutter.txt
Stutter Exceptions File	OCME_Fusion_Exceptions.csv
<i>Forward Stutter</i>	
Position Relative to Parent	1, 0
Inversely Proportional To	Expected Height of Stutter Peak
Stutter Regression File	OCME_Fusion_Forward Stutter.txt
Stutter Exceptions File	

8 Evaluation of the STRmix™ Analysis Diagnostics

- 8.1 The presence of a single sub-optimal diagnostic is not always an indication that rework is required. In some instances, a sub-optimal diagnostic(s) may be due to the nature of the sample (ex. low amounts of input DNA and/or stochastic effects) and not due to an issue with the STRmix™ run. Refer to the Troubleshooting Guide (Section 8) for further steps that may be taken to improve a sub-optimal diagnostic result.
- 8.2 For deconvolutions, verify that the following (**primary**) diagnostics conform to your qualitative expectations when compared to the electropherogram(s):
 - 8.2.1 The **mixture proportions** and **template amounts** assigned to the contributor(s).
 - 8.2.2 The **weights** assigned to the genotypes for each contributor listed in the **COMPONENT INTERPRETATION**.
 - 8.2.3 The **degradation** values and degradation plots (at the beginning of the report). The per contributor plot can be helpful for recognizing extreme degradation its potential influence over a single contributor's genotype weights across a profile.
 - 8.2.4 The **Locus Efficiencies (LSAE)**, following the variance charts.
- 8.3 For deconvolutions, evaluate the following (**secondary**) diagnostics for the run information listed in the **POST BURN-IN SUMMARY** at the beginning of the report.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 39 OF 43

- 8.3.1 **Total iterations:** If the total iterations exceeds 2.14 billion (2.14×10^9), this may lead to incorrect genotype weightings being assigned.
- 8.3.2 **Acceptance rate:** A very low acceptance rate (e.g. 1 in thousands to millions) may, in combination with the other diagnostics, indicate that the analysis needs to be run with additional accepts.
- 8.3.3 **Effective sample size (ESS):** A low ESS in relation to the total number of iterations suggests that the MCMC has not moved very far with each step or has had a low acceptance rate. A low ESS value (tens or hundreds) means that there is potential for a large difference in weights if the analysis is run again. For an ESS = 8, see the following troubleshooting guide.
- 8.3.4 **(Log)likelihood:** The larger this value, the better STRmix™ has been able to describe the observed data. A negative value suggests that STRmix™ has not been able to describe the data very well given the information it has been provided. A low or negative value for the log(likelihood) may indicate to users that the analysis requires additional scrutiny.
- 8.3.5 **Gelman-Rubin diagnostic:** If this value is above 1.2 then it is possible that the analysis has not converged, and the analysis requires additional scrutiny.
- 8.3.6 **LSAE variance:** The LSAE variance probability density chart in the report may be used in conjunction with the LSAE efficiency plot to identify unusual amplification within a profile (ex. extreme inhibition).
- 8.3.7 **Allele variance and stutter variance:** These variances should be compared to the modes and variance charts included in the reports. If the numbers are significantly elevated, the analysis may require additional scrutiny.
- 8.4 For LR comparisons, the overall **category of support** (inclusion, uninformative, exclusion) should conform to your qualitative expectations in comparison to the data.
- 8.4.1 Evaluate the **Per Locus Likelihood Ratio** table per locus and per sample, as well as the range of LR's **between population subgroups**; pay special attention to outliers and/or zero values.
- 8.5 For LR comparisons that result in support for an inclusion, check to ensure that the comparison sample falls in the appropriate contributor order. See report section titled **CONTRIBUTOR ORDER GIVING HIGHEST LR** at the beginning of the report.
- 8.5.1 If multiple comparison samples are positively associated with the same sample, results should be evaluated to ensure that they do not align with the same contributor. **See Section 9 for troubleshooting when more than one comparison aligns with the same contributor.**

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 40 OF 43

9 Troubleshooting Guide

- 9.1 The purpose of this guide is to address commonly seen scenarios which arise in casework. These guidelines are based on validation studies, literature references, and casework experience. However, not every situation can be covered by the Troubleshooting Guide. If a diagnostic issue arises that is not covered here, please discuss the issue with your supervisor, manager or the Technical Lead Team as needed.

Observations	Actions
<i>EVALUATION OF DECONVOLUTION DIAGNOSTICS</i>	
The mixture proportions or template amount do not reflect what is observed AND/OR the degradation does not reflect what is observed AND/OR the interpreted contributor genotypes do not conform to your qualitative expectations AND/OR the category of support for an LR comparison does not conform to your qualitative expectations	Re-evaluate the number of contributors; consider another STRmix™ analysis with one higher or one lower number of contributors. Consider amplifying a replicate if one has not already been done, with increased input amount when available and appropriate. Inhibition has occurred—microcon to clean and reamplify sample. Consider another STRmix™ analysis at a greater number of accepts (typically, 100,000 burn-in accepts and 500,000 total accepts per chain). Note: this requires approval by the Technical Leader (or his/her designee).
The mixture proportions and genotype weights do not conform to your qualitative expectations based on the electropherograms, and other troubleshooting options (including additional iterations) have been exhausted. Most common in mixtures with extreme ratios (ex. 98:2, 50:48:2,85:10:4:1) and/or many shared alleles between contributors)	Consider another STRmix analysis utilizing user informed Mx priors. This requires documented approval by the Technical Leader (or his/her designee).
The total iterations exceeds 2.14 billion (2.14 x 10 ⁹).	This could indicate the genotype weightings have been incorrectly assigned. Please contact the Technical Lead Team. An additional STRmix™ analysis may be required.
Effective sample size (ESS) = 8	Consider another STRmix™ analysis at a greater number of accepts (typically, 100,000 burn-in accepts and 500,000 total accepts per chain). Note: this requires approval by the Technical Leader (or their designee). You may want to use a computer with higher processing power.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 41 OF 43

A low or negative(log) likelihood	Re-evaluate number of contributors; consider another STRmix™ analysis with one higher or one lower number of contributors.
	Data has been removed that is allelic and/or stutter; data must be re-imported.
	Artifact peaks have been left labeled and must be removed.
Gelman-Rubin value is greater than 1.2	Consider another STRmix™ analysis at greater number of accepts (typically, 100,000 burn-in accepts and 500,000 total accepts per chain). This may sometimes reduce the GR to below 1.2. Note: this requires approval by the Technical Leader (or his/her designee).
Stutter and/or allele variance significantly elevated usually at or beyond the horizontal asymptote (may be in conjunction with low (log) likelihood)	Check to make sure no data has been omitted.
	Re-evaluate number of contributors; consider another STRmix™ analysis with one higher or one lower number of contributors.
	Consider amplifying a replicate if one has not already been done, with increased input amount when available and appropriate.
EVALUATION OF LR DIAGNOSTICS	
Differences in Per Locus LRs : large LR's (>1) obtained for each locus where the comparison profile is qualitatively included in the evidentiary profile, with one locus where the LR = 0	Data entry problem—check input files.
	Allele call not fully resolved at a given locus – ignore locus and perform analysis again. Note: discuss with supervisor as needed.
	Inhibition has occurred—microcon to clean and reamplify sample.
	Consider amplifying a replicate if one has not already been done, with increased input amount when available and appropriate.
	Re-evaluate number of contributors; consider another STRmix™ analysis with one higher or one lower number of contributors.
	Consider another STRmix™ analysis at greater number of accepts (typically, 100,000 burn-in accepts and 500,000 total accepts per chain). Note: this requires approval by the Technical Leader (or their designee).
Multiple elimination and/or comparison samples gave LRs supporting inclusion to an evidence sample and are aligning with the same contributor and/or there is an indication of relatedness .	Consult a supervisor and the Tech Lead Team to go over case specifics. Additional deconvolution and/or LR scenarios may need to be considered.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 42 OF 43

The 99% 1-sided HPD and Unified LR result(s) for one or more **population subgroup(s)** is several orders of magnitude lower than the other population subgroups in comparison to the point estimate LR, and the lowest LR of the four population subgroups does not conform to your qualitative expectations of the comparison.

Consult your supervisor or manager and the Technical Lead Team. A stratified likelihood ratio using census data may be calculated and reported; this must be approved by the Technical Leader (or their designee).

ERROR MESSAGES

Upon opening the STRmix software, an error that indicates “All connections to the floating license are currently in use”

Multiple licenses of STRmix™ are open on your computer. Close out the error message and determine if you have another active window of STRmix™ running.

Your previous use of the software was closed incorrectly. Re-start your computer, and if this does not resolve the issue, file an IT ticket to close your open STRmix™ license connections to the server.

All licenses are in use by other members of the laboratory.

Evidence profile data flag during pre-checks (as seen in step 3.16)

This is an example of a **missing stutter** flag. You will be notified if there’s an expected stutter peak that is not included based on the heights of the allelic peaks in the input file. See Section 3.16.1. Elevated stutter variance may occur in conjunction with this flag. The data should be evaluated to determine if the elevated variance makes intuitive sense with the peaks present in the profile. Consult the Tech Lead Team if further assistance is needed.

Another profile data issue that you may see is if there are peaks within the input file that are below the set analytical threshold. The original data should be re-evaluated in GeneMarker to determine if incorrect settings were applied.

Pre-Burnin failed:
Determine Acceptable Genotypes failed: Locus # (locus name) in the evidence cannot be explained given the parameters you have chosen.

Check your input value for NOC. Re-evaluate number of contributors; consider another STRmix™ analysis with one higher number of contributors.


Consider amplifying a replicate if one has not already been done, with increased input amount when available and appropriate.

(OutOfMemoryError) GC overhead limit exceeded

The run could not be completed with the computer power supplied. The run should be set up again using the Low Memory Mode setting and/or on a computer with more processing power.

FORENSIC BIOLOGY PROTOCOLS FOR FORENSIC STR ANALYSIS

STRmix v 2.7 Probabilistic Genotyping Standard Operating Instructions		
Status: Published		Document ID: 57029
DATE EFFECTIVE 05/04/2026	APPROVED BY Nuclear DNA Technical Leader	PAGE 43 OF 43

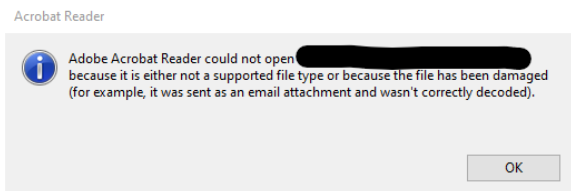

The results .pdf may be corrupted. Check that your run folder contains all other folders and files. The Reports folder may be missing the PNG and Excel files. You can regenerate your .pdf in STRmix by clicking on the Reports module from the main menu. Drag your full folder into the Previous Calculation field. Click the second Browse button that appears and add “_regen” to the end of your file name to not overwrite your original corrupted file. Include all files in the case folder. Add a case contact stating that a new .pdf was generated due to a corrupted file.

10 References:

- 10.1 STRmix™ v.2.7 Operation Manual and previous versions
- 10.2 STRmix™ v. 2.7 User’s Manual and previous versions
- 10.3 NYC OCME Internal Validation of STRmix™ v2.7 for Fusion 5C/3500xL (September 2021)
- 10.4 NYC OCME Internal Validation of STRmix™ v2.7 for Fusion 5C/3500xL Data - STRmix™ Parameters (August 2021)
- 10.5 NYC OCME Stutter Study for GeneMarker® HID 3.0.0 and STRmix™ Version 2.7- PowerPlex® Fusion Data run on 3500xl Genetic Analyzers (September 2021)
- 10.6 NYC OCME Internal Validation of STRmix™ V2.4 for Fusion (January 2017)
- 10.7 NYC OCME STRmix™ V2.4.08 Performance Check. (July 2018)
- 10.8 NYC OCME STRmix™ v2.7 Stutter Performance Check (April 2026)